

24- Well Ring Magnet Plate



Catalog No.: MSP24

Compatible With: 24-Well microplates
Round or Pyramid bottom

Used For: *Viscous samples,
DNA purification, and
nucleic acid purification*

Bead Type: *Any magnetic bead composition*

Material: *Solid aluminum design with hard-coat
finish for consistent results and maximum durability*

Automation Friendly: *ANSI/SBS compliant
footprint (127.75 x 85.50) with side clearance
for liquid handler grippers*

Working Volumes: 100 μ L - 10 mL

Magnets: *Dual Neodymium-Iron-Boron(NdFeB)
Ring magnet design*

Dimensions: W 127.75mm x L 85.50mm x H 15.9mm


Weight: 558 g / 1.23 lbs.

Country of Origin: USA

For research use only



Shown with 10 mL Whatman™ UniPlate™

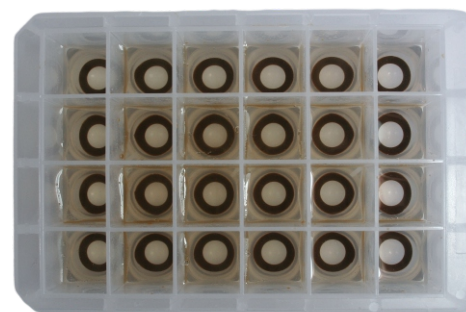

*Beads are separated into rings,
leaving a 10mm area at well bottom
for aspiration*



PERMAGEN®

email: info@permagenlabware.com
Tel. (978)914-4943

©2021 Permagen Labware LTD

Permagen Labware, LTD
300 Andover Street STE 169
Peabody, MA 01960
www.permagenlabware.com